

Insta Diagnostics Integration Guide

¹Table of Contents

Insta Diagnostics Integration Guide	1
Table of Contents	2
Insta Supported Diag Integration Messages	7
File Based Protocol	11
Send Orders (One Per Result label)	11
ORM	11
MSH Segment	13
PID Segment	14
PV1 Segment	15
IN1 Segment	15
IN2 Segment	16
DG1 Segment	17
ORC Segment	17
OBR Segment	18
Sending ORM (One Per Test)	18
ORM	18
MSH Segment	19
PID Segment	19
PV1 Segment	19
IN1 Segment (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)	19
IN2 Segment (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)	19
DG1 Segment (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)	19
ORC Segment (Below is the field which differ from above specification)	19
OBR Segment	19
Sending Cancellation	20
ORM	20
MSH Segment	20

PID Segment(This segment will be available based on the configuration 'send_pidpv1_for_cancel' in Insta)	21
PV1 Segment(This segment will be available based on the configuration 'send_pidpv1_for_cancel' in Insta)	21
IN1 Segment (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)	21
IN2 Segment (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)	21
DG1 Segment (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)	21
ORC Segment (Below is the field which differ from above specification)	21
OBR Segment	21
Sending Results	22
ORU	22
MSH Segment (Below is the field which differ from above specification)	22
PID Segment	22
PV1 Segment	22
ORC Segment (Below is the field which differ from above specification)	22
OBR Segment	22
OBX Segment	22
001	23
Sending Re-Conduction	24
ORM	24
MSH Segment	24
PID Segment	24
PV1 Segment	24
ORC Segment (Below are the fields which differ from above specification)	24
Sending Revert Sign-Off	25
ORM	25
MSH Segment	25
PID Segment	25
PV1 Segment	25
ORC Segment (Below are the fields which differ from above specification)	25
Sending Amended Results	25
ORM	25
MSH Segment	26
PID Segment	26
PV1 Segment	26
ORC Segment	26
NTE Segment	26

Sending Result As RP (Reference Pointer)	26
ORU	26
MSH Segment	27
PID Segment	27
PV1 Segment	27
ORC Segment	27
OBR Segment	27
OBX Segment (Below are the fields which differ from above specification)	27
Receiving Order	28
ORM	28
MSH Segment	28
PID Segment	28
PV1 Segment	28
ORC Segment	28
OBR Segment	28
Receiving Result (One Per Result label)	28
ORU	28
MSH Segment	28
PID Segment	28
PV1 Segment	28
ORC Segment	29
OBR Segment (Below are the fields which differ from above specification)	29
OBX Segment	29
Receiving Result As URL link	29
ORU	29
MSH Segment	30
PID Segment	30
PV1 Segment	30
ORC Segment	30
OBR Segment	30
OBX Segment (Below are the fields which differ from above specification)	30
Receiving Result As PDF/IMG	30
ORU	30
MSH Segment	31
PID Segment	31
PV1 Segment	31
ORC Segment	31
OBR Segment	31
OBX Segment (Below are the fields which differ from above specification)	31

Receiving Result (One Per Test)	31
MSH Segment	32
PID Segment	32
PV1 Segment	32
ORC Segment	32
OBR Segment (Below are the fields which differ from above specification)	32
OBX Segment	32
Receiving Re-conduction	33
ORM	33
MSH Segment	33
PID Segment	33
PV1 Segment	33
ORC Segment (Below are the fields which are differ from above specification)	33
OBR Segment	33
Receiving Results As Base64	33
ORU	33
MSH Segment	33
PID Segment	34
PV1 Segment	34
ORC Segment	34
OBR Segment	34
OBX Segment (Below are the fields which differ from above specification)	34
Receiving Revert Sign-off	34
ORM	34
MSH Segment	34
PID Segment	34
PV1 Segment	34
ORC Segment (Below are the fields which are differ from above specification)	35
OBR Segment	35
Receiving Amendment	35
ORM	35
MSH Segment	35
PID Segment	35
PV1 Segment	35
ORC Segment (Below are the fields which are differ from above specification)	35
OBR Segment	35
Receiving Cancellation	36
ORM	36

Socket Based Protocol	36
Send Order (One Per Result Label)	36
ORM	36
Sending ORM (One Per Test)	37
ORM	37
Sending Cancellation	37
ORM	37
Sending Addendum Sign -Off Report (Socket Based Protocol)	37
ORU	37
MSH Segment	38
PID Segment	38
PV1 Segment	38
ORC Segment (Below is the field which differ from above specification and Placer Group Number field is less)	38
OBR Segment	38
OBX Segment	38
Sending Commit Acknowledgement (Socket Based Protocol)	38
ACK	38
MSH Segment	38
MSA Segment	39
Sending Application Acknowledgement (Socket Based Protocol)	40
ACK	40
MSH Segment	40
MSA Segment	41
Receiving Addendum (Socket Based Protocol)	41
ORU	41
MSH Segment	42
PID Segment	42
PV1 Segment	42
ORC Segment (Below is the field which differ from above specification and Placer Group Number field is less)	42
OBR Segment	42
OBX Segment	42
Receiving Results (One Per Result Label, Socket Based Protocol)	43
ORU	43
MSH Segment	43
PID Segment	43
PV1 Segment	43
ORC Segment	43

OBR Segment	43
OBX Segment (Below are the supported formats in addition to the above specification)	43
Receiving Result (One Per Test)	43
ORU	43
Receiving Result As URL link	44
ORU	44
Receiving Results As Base64	44
ORU	44
Sending Sign-Off Report (Socket Based Protocol)	44
ORU	44
MSH Segment (Below are the fields which differ from above specification)	44
PID Segment (Below are fields which are differ from above specification)	44
PV1 Segment (Below are the fields which differ from above specification)	45
ORC Segment (Below is the field which differ from above specification and Placer Group Number field is less)	45
OBR Segment (Below are the fields which differ from above specification)	45
OBX Segment	45
DB Based Protocol	47
Sending DB Order	48
Order Entry	48
Sending DB Cancellation order	53
Receiving DB Results	53
Result Entry	53
Amendment and Re-conduction	54

Insta Supported Diag Integration Messages

The following types of messages are supported by Insta HMS:

Message	File Based Protocol		Socket Based Protocol		DB Based Protocol	
	Value	Template	Value	Template	Value	Template
Send Orders (one per Result Label)	Y (hl7 ORM message)	Not Applicable	Y (hl7 ORM message)	Not Applicable	Y Order (Non hl7)	Not Supported
Send Orders (one per Test)	Y (hl7 ORM message)	Y (hl7 ORM message)	Y (hl7 ORM message)	Y (hl7 ORM message)	Not Supported	Not Supported
Send Order Cancellation	Y (hl7 ORM message)	Y (hl7 ORM message)	Y (hl7 ORM message)	Y (hl7 ORM message)	Y Cancellation (Non hl7)	Not Supported
Send Re-Conducti on	Y (hl7 ORM message)	Y (hl7 ORM message)	Y (hl7 ORM message)	Y (hl7 ORM message)	Not Supported	Not Supported
Send Revert Signed-Off	Y (hl7 ORM message)	Y (hl7 ORM message)	Not Supported	Not Supported	Not Supported	Not Supported
Send Amendment	Y (hl7 ORM message)	Y (hl7 ORM message)	Not Supported	Not Supported	Not Supported	Not Supported
Send Results (Values)	Y (hl7 ORU message)	Y (hl7 ORU message)	Y (hl7 ORU message)	Y (hl7 ORU message)	Not Supported	Not Supported
Send Report as a file (pdf as a reference pointer)	Y (hl7 ORU message)	Y (hl7 ORU message)	Not Supported	Not Supported	Not Supported	Not Supported
Send Addendum (Text)	Not Supported	Not Supported	Not Supported	Y (hl7 ORU message)	Not Supported	Not Supported

Send Acknowledgment	Not Supported	Not Supported	Y (hl7 App Ack , hl7 Commit Ack)	Y (hl7 App Ack , hl7 Commit Ack)	Not Supported	Not Supported
Receive Order	Y (hl7 ORM message)	Y (hl7 ORM message)	Not Supported	Not Supported	Not Supported	Not Supported
Receive Cancellation	Y (hl7 ORM message)	Y (hl7 ORM message)	Not Supported	Not Supported	Not Supported	Not Supported
Receive Amendment	Y (hl7 ORM message)	Y (hl7 ORM message)	Not Supported	Not Supported	N Amendment (Non hl7)	Not Supported
Receive Re-Conduction	Y (hl7 ORM message)	Y (hl7 ORM message)	Not Supported	Not Supported	N Re-Conduction (Non hl7)	Not Supported
Receive Revert Signed-Off	Y (hl7 ORM message)	Y (hl7 ORM message)	Not Supported	Not Supported	Not Supported	Not Supported
Receive Result Values (one per Result label)	Y (hl7 ORU message)	Not Applicable	Y (hl7 ORU message)	Not Applicable	Y Result Receiving (Non hl7)	Not Supported
Receive Result (one per Test)	Y (hl7 ORU message)	Y (hl7 ORU message)	Y (hl7 ORU message)	Y (hl7 ORU message)	Not Supported	Not Supported
Receive Report as file (pdf/img stored on a separate folder(files))	Y (hl7 ORU message)	Y (hl7 ORU message)	Not supported	Not supported	Not Supported	Not Supported
Receive Report as ED (Encapsulate data)	Y (hl7_ORU_message)	Y (hl7_ORU_message)	Y (hl7_ORU_message)	Y (hl7_ORU_message)	Not Supported	Not Supported

Receive Report as URL link	Y (hl7 ORU message)	Y (hl7 ORU message)	Y (hl7 ORU message)	Y (hl7 ORU message)	Y Result Receiving (Non hl7)	Not Supported
----------------------------	---------------------------------------	---------------------------------------	---------------------------------------	---------------------------------------	--	---------------

Terminology

HL7 – Health Level 7 (<http://www.hl7.org/>)

Message – A collection of HL7 segments forming the basic unit of transmission.

Segment – A collection of fields describing one aspect of a complete HL7 message.

Field – One specific attribute of an HL7 Segment.

- File Based Protocol

1. Send Orders (One Per Result label)

ORM

The Segments/Groups in HL7 message include the following:

SEGMENT/GROUP	NAME	REQUIRED/OPTIONAL/REPEATABLE
MSH	Message header	Required
NTE	Notes and comments	Optional, Repeatable
Patient Group – Optional group (Required for new orders pertaining to a patient only)		
PID	Patient Identification	Required
PID1	Patient Demographics	Optional
NTE-1	Notes and Comments	Optional, Repeatable
Patient Visit Group – Optional group, part of Patient Group		
PV1	Patient Visit	Required
PV2	Patient Visit	Optional
Insurance Group – Optional and repeatable group, part of Patient Group		
IN1	Insurance	Optional
IN2	Insurance Additional Info	Optional
IN3	Insurance Additional Info Certification	Optional
Patient Group Continued		
GT1	Guarantor	Optional

AL1	Patient Allergy Information	Optional, Repeatable
Order Group - Repeatable Group		
ORC	Common Order	Required
Order Detail Group - Optional group, part of Order Group		
OBR	Observation Request	Required
NTE-2	Notes and comments	Optional, Repeatable
DG1	Diagnosis	Optional, Repeatable
Observation Group - Optional and Repeatable, part of Order Details Group		
OBX	Observation	Optional
NTE-3	Notes and comments	Optional, Repeatable
Order Group Continued		
CT1	Clinical trial identification	Optional
BLG	Billing	Optional

Sample ORM message

```
MSH|^~&|InstaHMS|HLM Tesing|||20161111154039||ORM^O01|6828411|P|2.4|||||^M
PID|1||MR000353^InstaHMS||Samplecollection^Cancellation^Bfr^^Captain||20081112|M|||Addres
s^^OTHER^002^^INDIA||99171878788^^Ruchi@gmail.com|||||||||^M
PV1|1|O|GENERAL PHYSICIAN^^GENERAL||||12017^ABDUL RAHIM.,M.D|||||DOC1275^Dr
ABDUL SAMEER||OP000487|||||||||^M
ORC|NW|2225617.4696.BLFPP000018^InstaHMS||494677|||||^M
```

OBR|1|2225617.4696.BLFPP000018^InstaHMS|2225617.4696.BLFPP000018|4696^List^InstaHMS^^|201611111532||||O||20161111153925122||DOC1275^Dr ABDUL SAMEER||BLFPP000018|10|||||||||||||20161111153200|^M

Field Separator	
Component Separator	^
Subcomponent Separator	&
Repetition Separator	~
Escape Character	\

MSH Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Field Separator	fieldSep	
2	Encoding Characters	encoding	^~\&
3	Sending Application (Sending/Receiving Application or Facility)	sendApp	InstaHMS
4	Sending Facility (Sending/Receiving Application or Facility)	sendFac	Diag OR Patient Registered Center Name (Incaset of multi center)
5	Receiving Application (Sending/Receiving Application or Facility)	recvApp	Nephrolife Care
6	Receiving Facility (Sending/Receiving Application or Facility)	recvFac	RBLRT

7	Date/Time Of Message	msgTimest	20161111154039
8	Security	security	
9	Message Type	msgType	ORM^O01
10	Message Control ID	controlld	12
11	Processing ID	proclld	P
12	Version ID	ver	2.4

PID Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Set ID - Patient ID	sid	1 (Always sends 1)
3	Patient ID (Internal ID) (includes Assigning Authority sub-component)	pidList	MR000353^InstaHMS (Unique id for a patient OR OldMno and Alternate)
5	Patient Name	name	Singh^Jitendar (patient's last name, patient name, middle name and salutation.)
7	Date/Time of Birth	dob	20081112
8	Sex	sex	M
11	Patient Address	addr	34, M Towers^^Bangalore^Karnata ka^^India
13	Phone Number - Home	phHome	9917187XXX^^^jitendar@gmail .com

			(Patient's phone and Email id)
--	--	--	--------------------------------

PV1 Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	SetID - PV1	sid	1
2	Patient Class	patientClass	O/I (Inpatient OR Outpatient)
3	Assigned Patient Location	assgndLoc	GENERAL PHYSICIAN^^GENERAL
7	Attending Doctor	attdDoctor	DOC002^Raghuram Bhat
8	Referring Doctor	refDoctor	12017^ABDUL RAHIM.,M.D
17	Admitting Doctor	admitDoctor	DOC1275^Dr ABDUL SAMEER
19	Visit Number	visitNum	OP000487
44	Admit Date/Time	admitTimest	20081112

IN1 Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
-----	---------------------	---------------------	---------

1	Set ID	sid	1
2	Insurance Plan ID	insPlanId	1233 (Numeric)
3	Insurance Company ID	insCold	ICM0009 (Alphanumeric)
4	Insurance Company Name	insCoName	Bajaj Allianz
5	Insurance Company Address	insCoAddr	Golden Heights,4th Floor,No.1/2, 59th C Cross Road, 4th M Block, 59A Cross Rd
7	Insurance Company Phone Number	insCoPhNum	078920 17448
12	Plan Effective Date	planEffectiveDate	2016111100000
13	Plan Expiration Date	planExpireDate	2017110100000
16	Name Of Insured	insuredName	Mohan
18	Insured's Date Of Birth	insuredDob	1980113000000
19	Insured's Address	insuredAddr	34, M Towers^^Bangalore^Karnataka^^India
43	Insured's Sex	insuredsAdministrativeSex	Male/Female

IN2 Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
25	Payor ID	payorId	MEM1233 (Member id)
61	Patient Member Number	patientMemberNum	12XXXXXX
63	Insured's Phone Number - Home	insuredsHomePhNum	080-XXXXXX

DG1 Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Payor ID	sid	1
3	Diagnosis Code	diagnosisCode	ARF
4	Diagnosis Description	diagnosisDescription	N17

ORC Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	NW
2	Placer Order Number	placerOrderNum	2225617.4696.BLFPP000018 ^InstaHMS

4	Placer Group Number	placerGrpNum	494677 (Insta Common Orderid)
---	---------------------	--------------	-------------------------------

OBR Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Set ID - OBR	sid	1
2	Placer Order Number	placerOrderNum	2225617.4696.BLFPP000018^InstaHMS
3	Filler Order Number	fillerOrderNum	2225617.4696.BLFPP000018
4	Universal Service ID	serviceId	901^LDL Cholesterol^InstaHMS
6	Requested Date/time	reqTime	201611111532
11	Specimen Action Code	specimenAction	O
13	Relevant Clinical Info	clinicalInfo	Data captured in Additional Order details section as user notes from the order screen.(Not supported for Internal lab tests)
14	Specimen Received Date/Time	specimenTime	20161111153925122
16	Ordering Provider	ordProvider	DOC1275^Dr ABDUL SAMEER
18	Placer field 1	placer1	BLFPP000018 (sample number)
19	Placer field 2	placer2	10 (sample type id)
36	Scheduled Date/Time	scheduledDtTime	20161111153200

2. Sending ORM (One Per Test)

ORM

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[IN1 Segment](#) (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)

[IN2 Segment](#) (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)

[DG1 Segment](#) (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)

[ORC Segment](#) (Below is the field which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
2	Placer Order Number	placerOrderNum	4696.^ V / T (Value Or Template)

OBR Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Set ID - OBR	sid	1
2	Placer Order Number	placerOrderNum	4696^ V / T
3	Filler Order Number	fillerOrderNum	4696

4	Universal Service ID	serviceId	Universal Service ID
6	Requested Date/time	reqTime	201611111532
11	Specimen Action Code	specimenAction	O
13	Relevant Clinical Info	clinicalInfo	Data captured in Additional Order details section as user notes from the order screen.(Not supported for Internal lab tests)
14	Specimen Received Date/Time	specimenTime	20161111153925122
16	Ordering Provider	ordProvider	DOC1275^Dr ABDUL SAMEER
36	Scheduled Date/Time	scheduledDtTime	20161111153200

3. Sending Cancellation

ORM

[MSH Segment](#)

[PID Segment](#)(This segment will be available based on the configuration 'send_pidpv1_for_cancel' in Insta)

[PV1 Segment](#)(This segment will be available based on the configuration 'send_pidpv1_for_cancel' in Insta)

[IN1 Segment](#) (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)

[IN2 Segment](#) (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)

[DG1 Segment](#) (This segment will be available based on the configuration 'Send Insurance Segments' in Insta)

[ORC Segment](#) (Below is the field which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	CA

[OBR Segment](#)

4. Sending Results

ORU

[MSH Segment](#) (Below is the field which differ from above specification)

9	R	Message Type	msgType	ORU^R01
---	---	--------------	---------	---------

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#) (Below is the field which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	RE

[OBR Segment](#)

OBX Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Set ID - OBX	sid	1

2	Value Type	valueType	TX / ST / RP / NM / NUM
3	Observation Identifier	obsId	CR002&BODY^RAD EX CH 1 VIEW FRONT
4	Observation Sub-ID	obsSubId	1
5	Observation Value	obsValue	7.0
6	Units	units	x10E3/uL
7	References Range	refRanges	3.4-10.8
8	Abnormal Flags	abnormalFlags	N / L / H / LL / HH
9	Probability	probability	
10	Nature of Abnormal Test	natureAbnTest	
11	Observ Result Status	obsResultSts	D, C OR F 001
12	Date Last Obs Normal Values	dateLastNormal	
13	User Defined Access Checks	accessChecks	
14	Date/Time of the Observation	obsTime	20151006181139
15	Producer's ID	producerId	Mediinterface
16	Responsible Observer	responsibleObserver	
17	Observation Method	obsMethod	

VALUE	LABEL
D	Deletes the OBX record
C	Record coming over is a correction and thus replaces a final result
F	Final result

5. Sending Re-Conduction

ORM

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#) (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	SC (Status changed)
5	Order Status	orderSts	RP (Order has been replaced)

6. Sending Revert Sign-Off

ORM

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#) (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	SC (Status changed)
5	Order Status	orderSts	IP (In process, Unspecified)

7. Sending Amended Results

ORM

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#)

NTE Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Set ID-NTE	sid	1
3	Comments	comment	Amendment reason
4	Comment Type	commentType	1R (Primary Reason)

8. Sending Result As RP (Reference Pointer)

ORU

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#)

[OBR Segment](#)

[OBX Segment](#) (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
2	Value Type	valueType	RP
3	Observation Identifier	obsId	24 Hr Urine,Echocardiography (Test names which involved in a report with comma separated)
5	Observation Value	obsValue	OP00012_0001.pdf (visitID_reportID)

9. Receiving Order

ORM

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#)

[OBR Segment](#)

10. Receiving Result (One Per Result label)

ORU

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

ORC Segment

OBR Segment (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
20	Filler Field 1	filler1	Supported any data type.In Insta it will be stored as a custom field, it has no value in Application,It is view purpose only, like in prints.
22	Results Rpt/Status Chng - Date/Time	resultTs	20161111154039

OBX Segment

11. Receiving Result As URL link

ORU

```
MSH|^~\&|MEDALL|DIAG|InstaHMS|DIAG|20160210125012||ORU^R01|||||
PID|||WLG22936^MEDALLHEALTHCARE||^VIVEK||19900101|F||^|9916624493|||||
PV1|1|O|||||9703^DR. PRAMOD.P||
ORC|RE|^MEDALLHEALTHCARE|||||
OBR|1||^MedallHealthCare|||||
OBX||RP|DGC0214^CT
Brain||http://www.hostname.com/report/result_report.pdf^^PDF^PDF||||F||20160702124717||^Dr
. KAMLESH|||^M
```

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#)

[OBR Segment](#)

[OBX Segment](#) (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
2	Value Type	valueType	RP
5	Observation Value	obsValue	http://hostname.com/result/ result_report.pdf

12. Receiving Result As PDF/IMG

ORU

```
MSH|^~\&|MEDALL|DIAG|InstaHMS|DIAG|20160210125012||ORU^R01|||||
PID|||WLG22936^MEDALLHEALTHCARE||^VIVEK||19900101|F||^|9916624493|||||
PV1|1|O|||||9703^DR. PRAMOD.P||
ORC|RE|^MEDALLHEALTHCARE|||||
OBR|1||^MedallHealthCare|||||
OBX|1|RP|^X-RAY CHEST
AP^MEDALLHEALTHCARE|1|1.2.392.200036.9125.9.0.253122479.2198387760.1537325039.jpeg
^RAD^Application^application/jpeg|SYS_RR|||||201602101237|
```

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#)

[OBR Segment](#)

[OBX Segment](#) (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
2	Value Type	valueType	RP
5	Observation Value	obsValue	1.2.392.200036.9125.9.0.2 53122479.2198387760.153 7325039.jpeg

13. Receiving Result (One Per Test)

ORU (With Multiple OBX's)

Sample message :

```
MSH|^~\&|InstaHMS|BRS|Nephrolife Care^2|PACS|20180221143829||ORU^R01|13894881|P|2.4
PID|1||MR012807^InstaHMS||-^Sadhana^^^Miss||19910101|F|||^^BANGALORE^017^^91||45354
PV1|1|O|General^^GENERAL|||DOC0189^Dr.Farha DeebaKhan|||DOC0232^ASMITHA||OP121045
ORC|NW|42586.2170^V||149656
OBR|1|42586.2170^V|42586.2170|2170^24 HR URINE PROTEIN/
ALBUMIN^InstaHMS||201802211402|||O|||DOC0232^ASMITHA||01|||01|||||||||20180221140200
OBX|1|NUM|6110060^WBC||6.3|x10E3/uL^x10E3/uL|3.4-10.8^3.4^10.8|LL|L|?
Unknown|F|||20171116160738||A.BRAVONRL^Bravo^Ana^^^^^^Personnel^^^External
Identifier|^^^CD:27807710
```

OBX|2|NUM|6110219^RBC||4.8|x10E6/uL^x10E6/uL|4.1-5.8^4.1^5.8|HH|H|?
 Unknown|F|||20171116160738||A.BRAVONRL^Bravo^Ana^^^^^^^Personnel^^^External
 Identifier|^^^CD:27807710
 OBX|3|NUM|6110061^Hemoglobin..||14.3|g/dL^g/dL|12.6-17.7^12.6^17.7|||?
 Unknown|F|||20171116160738||A.BRAVONRL^Bravo^Ana^^^^^^^Personnel^^^External
 Identifier|^^^CD:27807710
 OBX|4|NUM|43552^Hematocrit.||41.0|%^%%|37.5-51.0^37.5^51.0|||?
 Unknown|F|||20171116160738||A.BRAVONRL^Bravo^Ana^^^^^^^Personnel^^^External
 Identifier|^^^CD:27807710

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#)

[OBR Segment](#) (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
22	Results Rpt/Status Chng - Date/Time	resultTs	20161111154039

[OBX Segment](#)

14. Receiving Re-conduction

ORM

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#) (Below are the fields which are differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	SC (Status changed)
5	Order Status	orderSts	RP (Order has been replaced)

[OBR Segment](#)

15. Receiving Results As Base64

ORU

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#)

[OBR Segment](#)

[OBX Segment](#) (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
2	Value Type	valueType	ED
5	Observation Value	obsValue	^^PDF^BASE64^JVBERi0xLj QKJeLjz9MKMSA ...

16. Receiving Revert Sign-off

ORM

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#) (Below are the fields which are differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	SC (Status changed)
5	Order Status	orderSts	IP (In process, Unspecified)

[OBR Segment](#)

17. Receiving Amendment

ORM

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#) (Below are the fields which are differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	SC (Status changed)
5	Order Status	orderSts	IP (In process, Unspecified)

[OBR Segment](#)

18. Receiving Cancellation

[ORM](#)

- Socket Based Protocol

1. Send Order (One Per Result Label)

[ORM](#)

2. Sending ORM (One Per Test)

[ORM](#)

3. Sending Cancellation

[ORM](#)

4. Sending Addendum Sign -Off Report (Socket Based Protocol)

ORU

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#) (Below is the field which differ from above specification and Placer Group Number field is less)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	SC (Status Changed)

OBR Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
25	Result Status	resultSts	A (Some, but not all, results available)

[OBX Segment](#)

5. Sending Commit Acknowledgement (Socket Based Protocol)

ACK

MSH Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
-----	---------------------	---------------------	---------

1	Field Separator	fieldSep	
2	Encoding Characters	encoding	^~\&
3	Sending Application (Sending/Receiving Application or Facility)	sendApp	Insta
4	Sending Facility (Sending/Receiving Application or Facility)	sendFac	Diag
5	Receiving Application (Sending/Receiving Application or Facility)	recvApp	MWL
6	Receiving Facility (Sending/Receiving Application or Facility)	recvFac	MWL
7	Date/Time Of Message	msgTimest	20161111154039
8	Security	security	
9	Message Type	msgType	ACK
10	Message Control ID	controlId	12
11	Processing ID	procId	T
12	Version ID	ver	2.4
15	Accept Acknowledgement Type	ackType	NE (Never)
16	Application Acknowledgement Type	appAckType	NE (Never)

MSA Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Acknowledgement Code	code	CE / CA (Enhanced mode: Accept acknowledgment: Commit Accept / Enhanced mode: Accept acknowledgment: Commit

			Error)
2	Message Control ID	controlId	Control id from MSH
3	Text Message	msgText	Unparseable message / Unknown message type OR Empty message
4	Expected Sequence Number	expectedSequence	
5	Delayed Acknowledgement Type	delayedAckType	
6	Error Condition	errorCondition	

6. Sending Application Acknowledgement (Socket Based Protocol)

ACK

[MSH Segment](#)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
-----	---------------------	---------------------	---------

9	Message Type	msgType	ACK^R01
15	Accept Acknowledgement Type	ackType	AL (Always)

MSA Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Acknowledgement Code	code	AE / AA (Original mode: Application Accept - Enhanced mode: Application acknowledgment: Accept / Original mode: Application Error - Enhanced mode: Application acknowledgment Error)
2	Message Control ID	controlId	Control id from MSH
3	Text Message	msgText	Unparseable message / Unknown message type OR Empty message
4	Expected Sequence Number	expectedSequence	
5	Delayed Acknowledgement Type	delayedAckType	
6	Error Condition	errorCondition	0

7. Receiving Addendum (Socket Based Protocol)

ORU

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#) (Below is the field which differ from above specification and Placer Group Number field is less)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	SC (Status Changed)

OBR Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
22	Results Rpt/Status Chng - Date/Time	resultTs	20161111154039
25	Result Status	resultSts	A (Some, but not all, results available)

[OBX Segment](#)

8. Receiving Results (One Per Result Label, Socket Based Protocol)

ORU

[MSH Segment](#)

[PID Segment](#)

[PV1 Segment](#)

[ORC Segment](#)

[OBR Segment](#)

[OBX Segment](#) (Below are the supported formats in addition to the above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
2	Value Type	valueType	NM (Numeric)

9. Receiving Result (One Per Test)

[ORU](#)

10. Receiving Result As URL link

[ORU](#)

11. Receiving Results As Base64

[ORU](#)

12. Sending Sign-Off Report (Socket Based Protocol)

ORU

[MSH Segment](#) (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
4	Sending Facility (Sending/Receiving Application or Facility)	sendFac	Diag OR Center Code and Patient Registered Center Name (Incaset of multi center)
9	Message Type	msgType	ORU^R01

[PID Segment](#) (Below are fields which are differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
5	Patient Name	name	Singh^Jitendar (patient's last name, patient name and salutation.)

13	Phone Number - Home	phHome	9917187XXX^^^jitendar@gmail.com (Patient's phone)
----	---------------------	--------	--

[PV1 Segment](#) (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
2	Patient Class	patientClass	O/I(Inpatient, Incoming Patient OR Outpatient)

[ORC Segment](#) (Below is the field which differ from above specification and Placer Group Number field is less)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Order Control	ctrl	RE (Observations/ Performed Service to follow)


[OBR Segment](#) (Below are the fields which differ from above specification)

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
25	Result Status	resultSts	F (Final Result)

OBX Segment

SEQ	ELEMENT DESCRIPTION	INTERNAL IDENTIFIER	EXAMPLE
1	Set ID - OBX	sid	1
2	Value Type	valueType	TX
3	Observation Identifier	obsId	(HL7Exportcode / testid and Testname)
5	Observation Value	obsValue	Report content

- DB Based Protocol

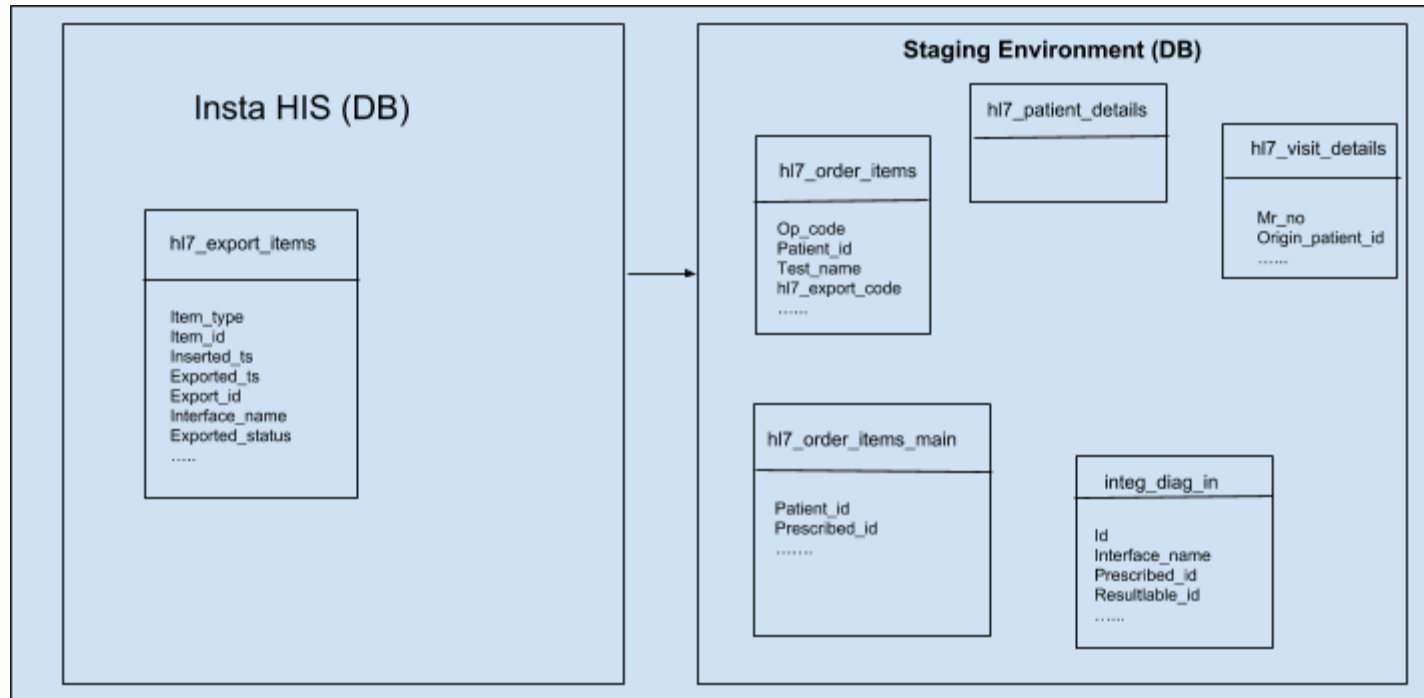

1. Sending DB Order

Order Entry

Staging table :hl7_order_items

DB Column	EXAMPLE
item_type	TEST,TESTMODIFIED,TESTRESULT
item_id	45629(unique identification number in insta for each order)
inserted_ts	2017-10-03 15:47:46.845976
exported_ts	2017-10-03 17:40:46.642376
export_id	32435(Unique export sequence)
interface_name	medi_interface
export_status	N-not yet exported, F-failed to export, S-export successful
export_msg_id	14322(message sequence)
export_failure_msg	Failure message
bill_paid	N-bill is not paid or sample not collected, Y- bill is paid
op_code	N-order to read, C- order cancelled indicating not to read
patient_id	OP120725(visit id)
test_name	COMPLETE BLOOD COUNTS (CBC)
hl7_export_code	
test_id	DGC0783
pres_date	2017-01-03 15:20:00
common_order_id	43555(order number)
scheduled_time_start	2013-09-13 15:24
start_date_time	2013-09-13 15:30:10
scheduled_time_end	2013-09-13 17:24

equipment_code	
sample_sno	SE10911
sample_date	2017-01-03 15:20:00(sample collected date)
mr_no	MR013140
oldmrno	PR14065
last_name	xyz(patient last name)
patient_name	abc(patient full name)
salutation	Mr, Ms
patient_gender	M-Male, F-Female
expected_dob	19830404 (expected date of birth calculated by age)
patient_address	BANGALORE
cityname	BANGALORE
city_code	088
statename	KARNATAKA
state_code	017
country_name	India
country_code	91
patient_phone	7800938764
reg_date	2016-03-30
doctor	DOC0189(unique doctor id for each doctor from insta)
doctor_name	DR.SUSHMA RANI
bed_type	GENERAL1,GENERAL2
ward_name	GENERAL
reference_docto_id	DOC0067 (refferal doctor id)
diagnosis_code	N17.8 (icd_code)
diagnosis	Hypercholesterolemia, metabolic syndrome

visit_type	OP-out patient, IP- in patient, OSP- outside patient
refdoctorname	DR.NAMITHA (refferal doctor name)
conduction_format	T-Template based test, V-Value based test
resultlabel_id	233(result id)
resultlabel	Hemoglobin
center_id	5(unique center id from insta)
center_name	Koramangala region
conducting_doctor_id	DOC0189(conducting doctor id)
conducting_doctor_name	DR.SUSHMA RANI (conducting doctor name)
prescription_doctor_id	DOC0189(prescription doctor id)
prescription_doctor_name	DR.SUSHMA RANI (test prescribed doctor name)
prescription_department	Microbiology (prescribed doctor department)
out_house_sample	t- true(out_house test),f- false(not out-house test)
middle_name	Mn (patient middle name)

Staging table: HI7_order_items_main

DB Column	EXAMPLE
patient_id	IN000250
prescribed_id	67758(item_id)
inserted_ts	2017-10-03 15:47:46.845976
exported_ts	2017-10-03 17:40:46.642376
export_id	32435(Unique export sequence)
package_name	Body check up (package name)
tpa_name	BAJAJ ALLIANZ
plan_name	Bajaj plan
employee_id	P3244

order_remarks	Order details
additional_test_info	Clinical information
orig_center_id	7
priority	R-Regular, S-Stat, U-Urgent
sample_type	Blood,serum,urine etc..
email_id	abc@gmail.com
origin_center_name	Koramangla lab center
origin_patient_id	OP121650
out_house_name	t- true(out_house test),f- false(not out-house test)

Insurance modified details available in **hl7_visit_details** staging table (Visit wise)

DB Column	EXAMPLE
mr_no	MR013140
origin_patient_id	OP121650
inserted_ts	2017-10-03 15:47:46.845976
exported_ts	2017-10-03 17:40:46.642376
tpa_name	BAJAJ ALLIANZ
plan_name	Bajaj plan
employee_id	P3244
patient_details_status	'VM'- Visit modifications

Patient demographic changes available in **hl7_patient_details** staging table (mr_no wise)

DB Column	EXAMPLE
hl7_patient_details_id	233(sequence no)
mr_no	MR013140
oldmrno	PR14065
last_name	xyz(patient last name)
patient_name	abc(patient full name)
salutation	Mr, Ms
patient_gender	M-Male, F-Female
expected_dob	19830404 (expected date of birth calculated by age)
patient_address	BANGALORE
cityname	BANGALORE
city_code	088
statename	KARNATAKA
state_code	017
country_name	India
country_code	91
patient_phone	7800938764
inserted_ts	2017-10-03 15:47:46.845976 (Timestamp when order got generated in Insta)
exported_ts	2017-10-03 17:40:46.642376 (Timestamp when Insta exports to staging table)
email_id	abc@gmail.com
patient_details_status	'DM' (Indicates patient demography modified)

2. Sending DB Cancellation order

[Sending DB Order](#)(Below is the field which differ from above specification of hl7_order_items staging table to initiate order cancellation from insta)

DB Column	EXAMPLE
op_code	C (order cancellation from insta)

3. Receiving DB Results

Result Entry

Staging table to read results in insta: **integ_diag_in**

DB Column	EXAMPLE
id	3243(Sequence number)
interface_name	medi_interface
prescribed_id	34567(Unique item id for each order)
resultlabel_id	34(Result Label id number)
sample_no	SE10911 (sample identification number)
report_value	2.3(result label value)
status	N(normal)- result to read, RR(report ready)- external result to read, SC(status changed)- amendment/reconduction request, S- successfully read, F-Failed to read.
comments	Result entry comments
inserted_ts	2017-10-03 15:47:46.845976(result entry time)
read_ts	2017-10-03 15:50:46.845976(results read time)
export_failure_msg	Failure message

value_type	ST-value based test, TX- template based test
pat_id	Insta patient_id (visit id)

4. Amendment and Re-conduction

[Receiving DB Results](#)(Below is the field which differs from above specifications of integ_diag_in table to initiate amendment/reconduction)

DB Column	EXAMPLE
status	SC(status changed) - insta supports amendment/Re-conduction.